Inbetriebnahme / Laden
Mise en service / Charge
Messa in funzione / Caricamento
Initial operation / charging التشغيل / الشحن

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Read this manual carefully before first use and keep it for future reference. Use the device exclusively for cleaning teeth. The safety instructions must be observed at all times during use, storage and disposal. Electrical specification: see bottom of charging device.

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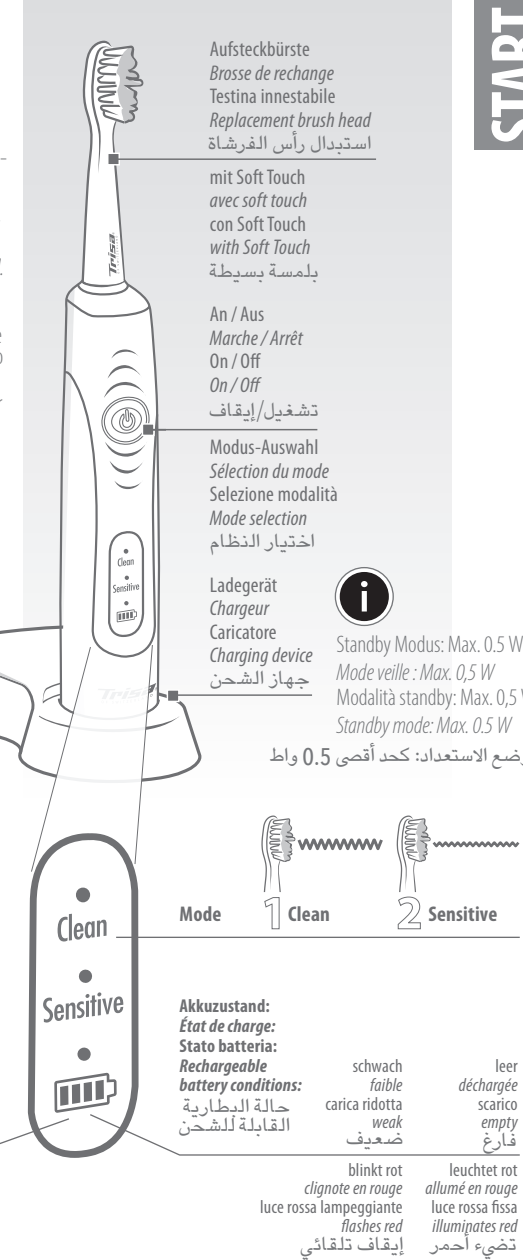
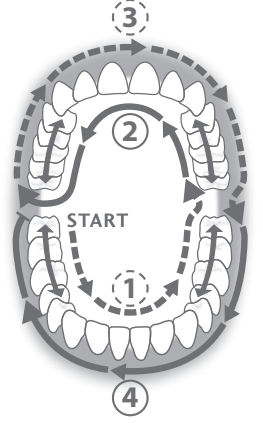
A brief stop provides a reminder after every 30 seconds to change over to the next tooth area. Thereby, you are certain to have brushed each section for a sufficient period of time. After 2 minutes, a longer stop signalizes the end of tooth brushing.

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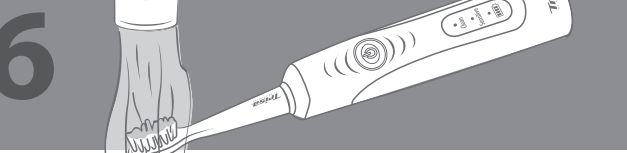
• Replace brush head every 3 months.
• For hygiene / functionality, rinse brush head and switched-on appliance under running water after each use.

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Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.

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Check supply cord and charging device regularly for cracks or other damages. Never touch with sharp/hot objects. The supply cord cannot be replaced. If the cord is damaged the appliance should be scrapped. Never operate damaged charging device or toothbrush. Improper repairs or modifications to the device may result in electric shock, injury or fire.

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- **WARNING:** For the purposes of recharging the battery, only use the detachable supply unit 4680, 4656 or 4663 provided with this appliance.

Battery Exchange Only have the battery replaced by qualified specialists. Use Li-ion battery type 14500, 3.7V, min. 750mAh with connecting wires. Use a coin to open the cover by turning it a quarter turn clockwise. Carefully push out the inner part by pressing on the steel axle. Check the battery polarity to ensure correct installation of the new battery. Use a soldering iron to disconnect the connection wires of the old battery from the electronics. Carefully remove the battery from the holder. Insert the new battery into the holder and solder the connection wires with the correct polarity. Position the seal precisely on the inner part. Carefully insert the inner part so that the seal is flush with the housing outlet. Replace the cover and close it again.

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The manufacturer will provide you with a two-year guarantee for this appliance from the date of purchase (purchase receipt). The guarantee service shall comprise a repair or replacement of an appliance with material or manufacturing defects. In order to maintain the guarantee, TRISA original brushes shall exclusively be used. The appearance of wear and tear within the context of normal use shall not be subject to the guarantee. The guarantee shall not be subject to the guarantee. The guarantee shall not be subject to the guarantee.

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posed of with household waste. Make use of the return possibility in the designated shops or collection centres.

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tee shall lapse in case of third-party intervention or with changes to the original condition, as well as commercial use. The guarantee performance requires that the defective appliance be submitted with the purchase receipt.

Hinweis für die Entsorgungsstelle Remarque pour le centre de récupération Indicazioni per il centro di raccolta per lo smaltimento Information for the disposal centre

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This toothbrush contains a lithium ion rechargeable battery. For environmental reasons, it must be removed and disposed of in an environmentally compatible manner at the end of its service life.

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- TRISA AG Kantonsstrasse 31 | CH-6234 Triengen Switzerland info@trisa.ch | 041 935 35 35 www.trisa.ch